


# Clear-IT™

## Sample Marvelously with Elcam's new Closed Blood Sampling System

The Clear-IT™ Closed Blood Sampling System is an advanced tool designed by Elcam in response to two common problems associated with blood sampling: ICU Anemia due to blood waste associated with blood sampling, and infections due to exposure of patient's blood system to the atmosphere.

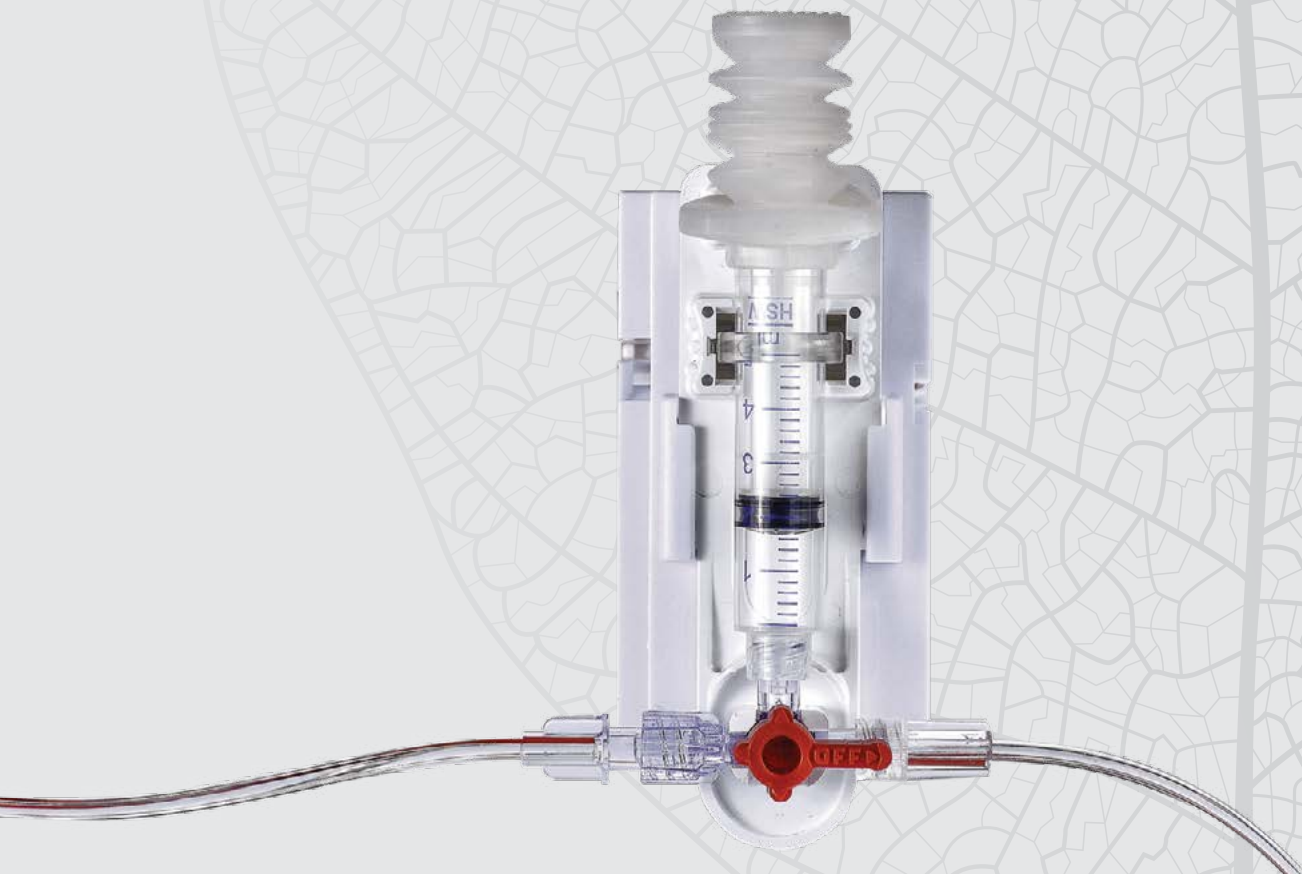

Patient Monitoring




**Clear-IT™** Closed Blood Sampling System is designed for simplicity, safety and ease of use.

The **Clear-IT™** combines an **In-Line Reservoir** for waste collection and a **stopcock with Marvelous™ handle** that enables flushing the line and syringe tip with minimal residual volume.

The **In-Line Reservoir** eliminates the discarding of waste taken through the set prior to the blood sampling. The entire clearing volume is reinfused to the patient upon completion of the sampling process in a safe manner that minimizes contamination. The Reservoir is designed to be connected to a standard mount plate.

The **stopcock with the Marvelous™ proprietary handle design** and its circumferential channel, assures clearing of the line and the syringe's tip with minimal unflushable (residual) volume after blood reinfusion following blood sampling.

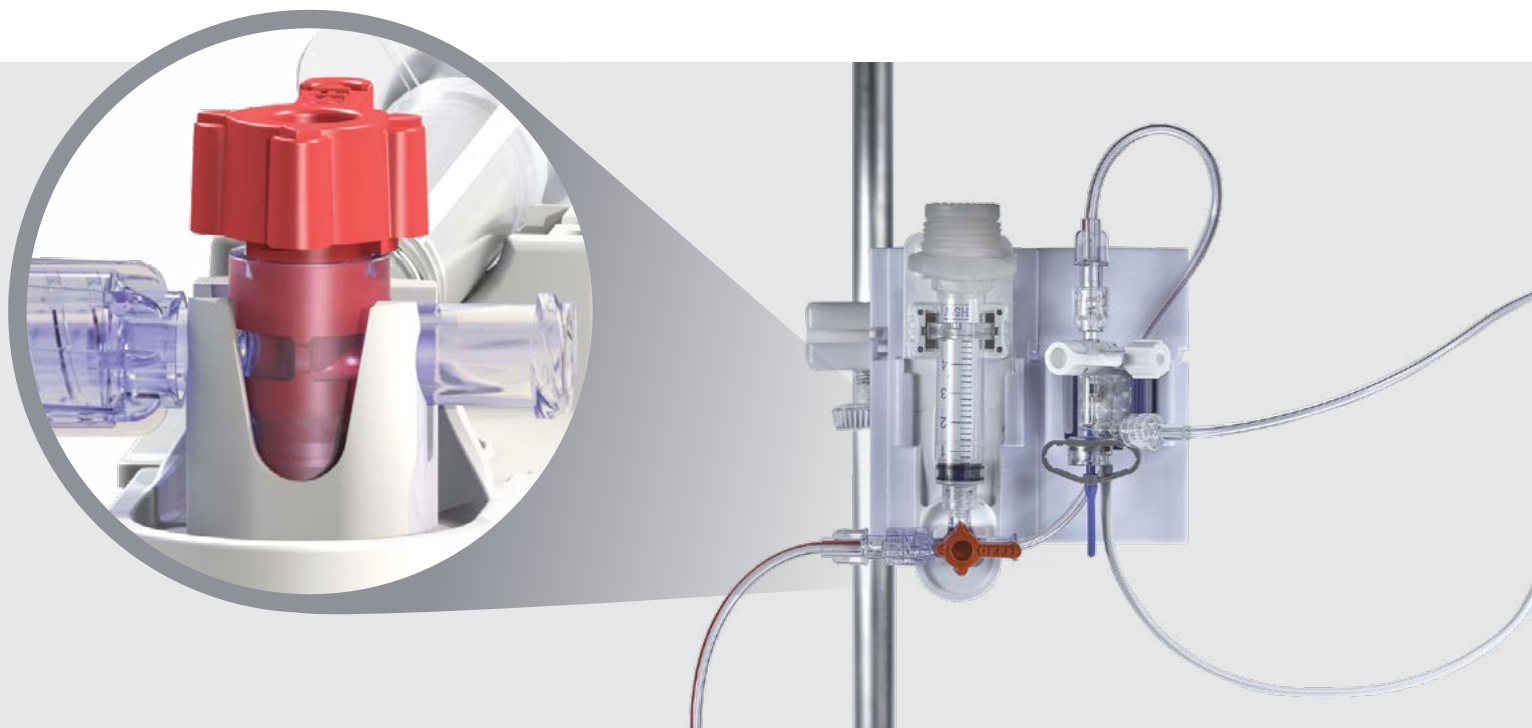

## Benefits

- Intuitive & Easy to use - Minimal training required
- Efficient - Reduces blood waste by reinfusing clearing volume back into the patient
- Safe - (\*) Closed in-line design protects patient's blood from exposure to bacteria  
(\* ) Continuous line and syringe clearing decreases opportunities for bacterial contamination

The Clear-IT™ is available only as an OEM product

**Elcam's Clear-IT™ Blood Sampling System is available as a stand-alone system. The system can be used in combination with Elcam's Sense- IT™ Integrated Disposable Pressure Transducer and with Elcam's Marvelous™ stopcock for minimal residual volume needless blood sampling.**



[www.elcam-medical.com](http://www.elcam-medical.com)

Head Office / Manufacturer  
Elcam Medical

BarAm 1386000,  
Israel  
Tel: 972 4 698 8120/1  
Fax: 972 4 698 0777  
sales@elcam.co.il

USA  
Elcam Medical Inc.

2 University Plaza, Suite 620  
Hackensack, NJ 07601  
Tel: +1 201 4571120  
Fax: +1 201 4571125  
info@elcam-medical.com

Europe  
Elcam Medical Europe

Via Emilia Romagna 15,  
41012 Carpi (Modena) Italy  
Tel: +39 059 697710  
Fax: +39 059 697705  
sales.eu@elcam-medical.com

China  
Elcam China Office

Room 3309, Jinhai Mansion  
No. 89 of Anshun Road,  
Shanghai, PRC. Zip: 200051  
Tel: +86 21 62829967  
Fax: +86 21 62829967  
johnliu02@126.com