


Elcam Stopcocks


Closed Swabbable Stopcock

A Closed Stopcock, integrated with a Swabbable Luer Activated Valve, allowing for maintenance of a fully closed system through the entire use of the stopcock, as well as eliminating needle stick injuries when using a stopcock.

Elcam's Closed Swabbable Stopcock is part of our E-Safe™ product family, specially developed in order to provide our customers with a range of solutions to minimize stopcock contaminations and enhancing patient safety.

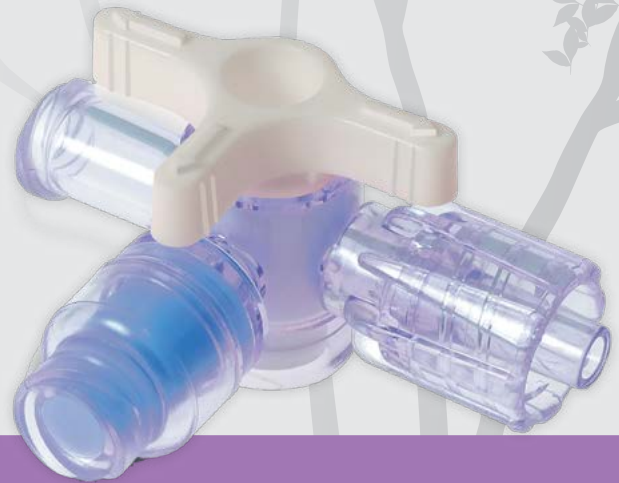


IV Therapy & Anesthesia

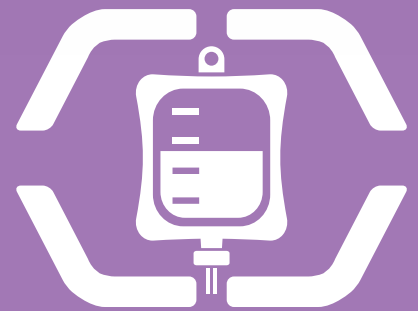

Benefits

Safety


- Normally Closed luer activated valve - Provides full barrier protection against pathogens entering through the stopcock ports
- Swabbable valve - allows effective disinfection of valve surface
- Needle free connection - prevents needle-stick injuries
- Low profile - minimizes dead space in comparison to capping with additional LAV

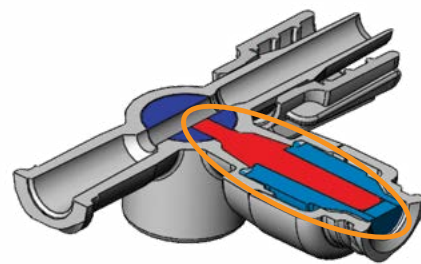

Residual volume in closed swabbable stopcock design

Ease of use

- Replaces the need to add a cap or Luer Activated Valve on a stopcock
- Eliminates the need for additional adaptors or canula

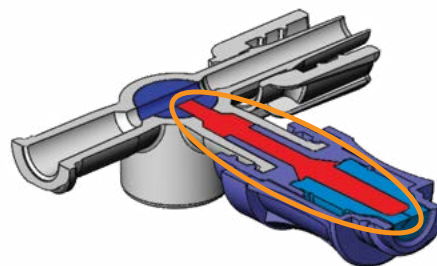

Residual volume in standard stopcock with LAV cap

Applications

- Suitable for all IV lines and anesthesia applications
- Multilumen Central Venous Catheters

Physical Properties

- Residual volume at the valve port - 0.12cc (High Flow)
- Minimum flow rate of 150 ml/min under 1 meter water column height
- The device is biocompatible according to ISO 10993-1
- Lipid Resistant
- Not made with latex
- Non-hemolytic

Configuration options

- 1 Valve or 2 valves
- Handles: One arm or 3 arm
- Different Handle colors
- High Flow or Standard Flow
- With or without Luer Lock
- Manifold

The Closed Swabbable Stopcock is CE approved for non-sterile products and 510(k) cleared



www.elcam-medical.com

Head Office / Manufacturer
Elcam Medical

BarAm 1386000,
Israel
Tel: 972 4 698 8120/1
Fax: 972 4 698 0777
sales@elcam.co.il

USA
Elcam Medical Inc.

2 University Plaza, Suite 620
Hackensack, NJ 07601
Tel: +1 201 4571120
Fax: +1 201 4571125
info@elcam-medical.com

Europe
Elcam Medical Europe

Via Emilia Romagna 15,
41012 Carpi (Modena) Italy
Tel: +39 059 697710
Fax: +39 059 697705
sales.eu@elcam-medical.com

China
Elcam China Office

Room 3309, Jinhai Mansion
No. 89 of Anshun Road,
Shanghai, PRC. Zip: 200051
Tel: +86 21 62829967
Fax: +86 21 62829967
johnliu02@126.com