

Elcam Stopcocks

HIGH FLOW STOPCOCKS

High Flow Path Ways

Designed mainly to provide higher flow rates or viscous fluid applications where required. Major applications are within rapid infusion, anesthesia and blood related procedures.



Your source of consistent quality and care








HIGH FLOW STOPCOCKS

High Flow Path Ways

Elcam's High Flow Stopcocks are available in 1,3 and 4 way configurations with various handle styles and luer locks to match different required applications. Our High Flow Stopcocks can be supplied in materials suitable for ETO and Gamma irradiation, and/or lipid based admixtures upon request.

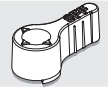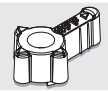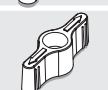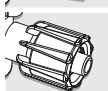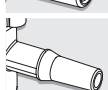

			1 way (90°/360°)	3 way (180°)	4 way (360°)
Handle Type		High Profile		<input type="checkbox"/>	<input type="checkbox"/>
		Low Profile		<input type="checkbox"/>	<input type="checkbox"/>
		GXL Handle		<input type="checkbox"/>	<input type="checkbox"/>
		T - Handle			<input type="checkbox"/>
		1 Arm HF	<input type="checkbox"/>		
		2 Arm High	<input type="checkbox"/>		
End Connections		Luer Lock	<input type="checkbox"/>		
		Short Luer Lock		<input type="checkbox"/>	<input type="checkbox"/>
		Long Luer Lock		<input type="checkbox"/>	<input type="checkbox"/>
		Fixed Luer Lock		<input type="checkbox"/>	<input type="checkbox"/>
		Luer Slip	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
		Tube Pocket	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

Our Stopcocks can be supplied with or without vented or sealed port covers.

At Elcam Stopcocks you can find the widest variety of "off the shelf" stopcocks tubing ports, for a variety of tubes, at all ports.

Body Material:

Polycarbonate or Polysulfone.

Raw material available: either clear (for EtO sterilization only) or tinted blue (for gamma).

Colors:

Subject to stopcock configuration.



Quality in every aspect:

Elcam's highly automated processes ensure production of consistently high quality products. We operate in strict compliance with the international quality standards. Elcam's High Flow Stopcocks are CE approved for non-sterile products and 510(k) cleared.



www.elcam-medical.com

Head Office / Manufacturer
Elcam Medical

BarAm 1386000,
Israel
Tel: 972 4 698 8120/1/2
Fax: 972 4 698 0777
sales@elcam.co.il

USA
Elcam Medical Inc.

2 University Plaza, Suite 620
Hackensack, NJ 07601
Tel: +1 201 4571120
Fax: +1 201 4571125
info@elcam-medical.com

Europe
Elcam Medical Europe

Via Emilia Romagna 15,
41012 Carpi (Modena) Italy
Tel: +39 059 697710
Fax: +39 059 697705
sales.eu@elcam-medical.com

China
Elcam China Office

Room 3309, Jinhai Mansion
No. 89 of Anshun Road,
Sanghai, PRC. Zip: 200051
Tel: +86 21 62829967
Fax: +86 21 62829967
johnliu02@126.com