# Elcam Stopcocks Closed Swabbable Stopcock

A Closed Stopcock, integrated with a Swabbable Luer Activated Valve, allowing for maintenance of a fully closed system through the entire use of the stopcock, as well as eliminating needle stick injuries when using a stopcock.

Elcam's Closed Swabbable Stopcock is part of our E-Safe™ product family, specially developed in order to provide our customers with a range of solutions to minimize stopcock contaminations and enhancing patient safety.

IV Therapy & Anesthesia





### Safety

- ➤ Normally Closed luer activated valve Provides full barrier protection against pathogens entering through the stopcock ports
- > Swabbable valve allows effective disinfection of valve surface
- ➤ Needle free connection prevents needle-stick injuries
- ➤ Low profile minimizes dead space in comparison to capping with additional LAV

### Ease of use


- Replaces the need to add a cap or Luer Activated Valve on a stopcock
- > Eliminates the need for additional adaptors or canula

# **Applications**

- > Suitable for all IV lines and anesthesia applications
- ➤ Multilumen Central Venous Catheters


### **Physical Properties**

- > Residual volume at the valve port 0.12cc (High Flow)
- ➤ Minimum flow rate of 150 ml/min under 1 meter water column height
- ➤ The device is biocompatible according to ISO 10993-1
- ➤ Lipid Resistant
- > Not made with latex
- ➤ Non-hemolytic

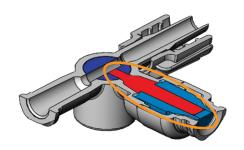

REV-8A

Residual volume in closed swabbable stopcock design

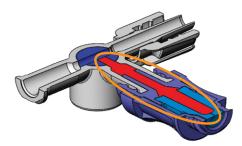

Residual volume in standard stopcock with LAV cap

## **Configuration options**

- ➤ 1 Valve or 2 valves
- ➤ Handles: One arm or 3 arm
- ➤ Different Handle colors
- ➤ High Flow or Standard Flow With or without Luer Lock
- **≻** Manifold

The Closed Swabbable Stopcock is CE approved for non-sterile products and 510(k) cleared



Head Office / Manufacturer Elcam Medical

BarAm 1386000, Israel Tel: 972 4 698 8120/1 Fax: 972 4 698 0777 sales@elcam.co.il USA Elcam Medical Inc.

1 University Plaza, Ste. 511 Hackensack, NJ 07601 Tel: +1 201 4571120 Fax: +1 201 4571125 info@elcam-medical.com Europe Elcam Medical Europe

Via Emilia Romagna 15, 41012 Carpi (Modena) Italy Tel: +39 059 697710 Fax: +39 059 697705 sales.eu@elcam-medical.com China Elcam China Office

Room 3309, Jinhai Mansion No. 89 of Anshun Road, Shanghai, PRC. Zip: 200051 Tel: +86 21 62829967 Fax: +86 21 62829967 johnliu02@126.com