




Dialysis Product Line

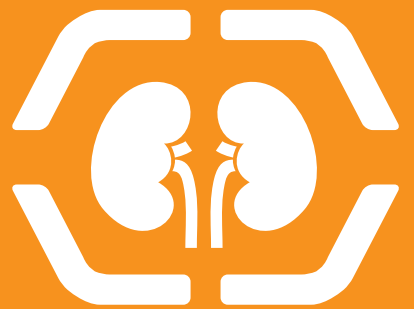

 [Elcam Devices](#)

 [Elcam Stopcocks](#)


 [Elcam Accessories](#)

Dialysis

As the Dialysis population continues to grow, treatment solutions whether In-center dialysis, home care dialysis and others, continue to evolve and improve through focusing on devices that provide the best results for these patients. Elcam Medical's Dialysis product line comprise a wide selection, providing our customers with solutions for their blood line products, whether customized or standard. Our devices and accessories are designed for easy and convenient operation and for promoting patient safety



Elcam Devices for Dialysis

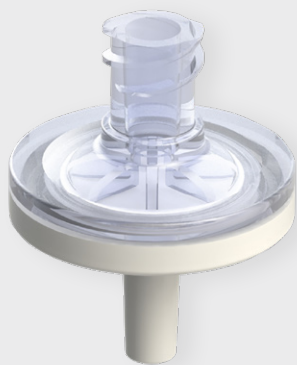

Transducer Protectors

Elcam's wide selection of transducer protectors grant the perfect barrier to contamination deriving from dialysis blood line applications.

We offer a comprehensive selection of transducer protectors designed for a simple and easy connection to the hemodialysis machine. Ranging from simple double tube fitting, to FLL/MLL or luer port and snap-off connector, as well as air vent with dedicated membrane. Our automated assembly technology performs a 100% in-process testing for tightness of the ultrasonic welding.



Drip Chambers

Drip Chambers are an essential part of the dialysis blood line set. They need to feature high quality and convenience, as well as enable customization according to different requirements.

Elcam's dialysis drip chambers incorporate our high quality and service standards. They are available in different diameters and tube ports, for pediatric and adult applications, including a large range of accessories like wall adaptors, fluid filters, lids with tube pocket and with male luer lock.

Elcam Accessories for Dialysis

Elcam's Accessories for Dialysis applications provide you with an extensive product offering designed to fit the everyday needs of dialysis fluid control applications. Our components can be customized according to you and your customers' needs, and can be supplied in a variety


of materials validated for EtO and gamma sterilization, not being made with DOP and DEHP. Available in different colors, tube fitment sizes and more.

Among our specialty solutions:



Safe 2 Rotator™ Luer Connector

Prevents accidental disconnection and tube kinking and provides secured connection along the entire treatment to reduce associated risks like blood loss, drug spills and air ingress. In compliance with ISO 80369-7.

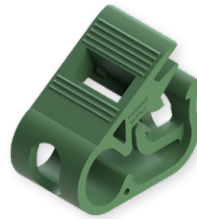

Non-Reopening Clamps

Designed for situations where it is imperative that the clamping device remains closed to avoid leaks or contamination, as in the case of dialysis solution bags.



Pressure Pillows


Soft reservoir allowing for reduction of pressure leaps in fluid circuits. The Elcam pillows are available in standard version and with non-return valve, with service line port.

Elcam offers many more accessories for dialysis applications:

Male and female Luer locks, connectors, injection sites, caps, clamps and more.

All these products are provided with our high-quality standards as well as Elcam's renowned excellent service and support, from delivery time to design changes, from biocompatibility to optimizing bonding processes.

Typical Arterial and Venous Blood Line Sets for Dialysis

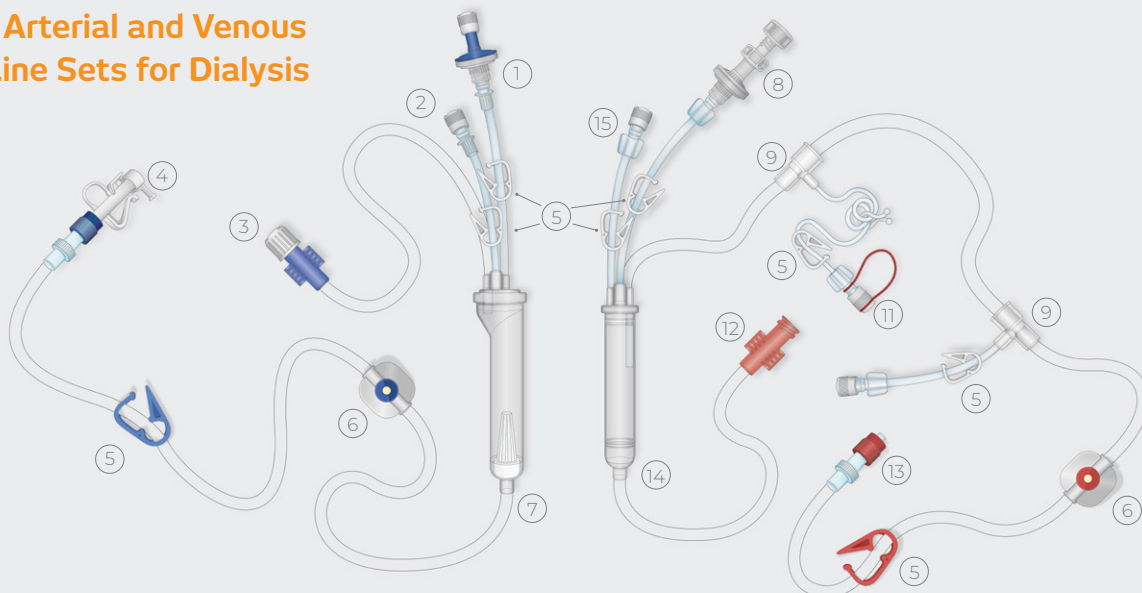

- | | | |
|---|--|---|
| 1. Transducer Protector | 6. Injection Sites | 11. FLL Connector plus Cap with fixing ring |
| 2. FLL Connector plus Caps | 7. Venous Drip Chamber with Filter | 12. Dialyzer Connector |
| 3. Dialyzer Connectors with Cap | 8. Transducer Protector with soft snap-off and Connector | 13. Arterial Patient Connector |
| 4. Recirculating Connector + Venous Patient Connector | 9. Pump Connectors | 14. Arterial Drip Chamber |
| 5. Clamps | 10. Heparine Clip | 15. FLL Connectors plus Caps |

**Head Office / Manufacturer
Elcam Medical**

BarAm 1386000,
Israel
Tel: 972 4 698 8120/1
Fax: 972 4 698 0777
sales@elcam.co.il

**USA
Elcam Medical Inc.**

1 University Plaza, Ste. 511
Hackensack, NJ 07601
Tel: +1 201 4571120
Fax: +1 201 4571125
info@elcam-medical.com

**Europe
Elcam Medical Europe**

Via Emilia Romagna 15,
41012 Carpi (Modena) Italy
Tel: +39 059 697710
Fax: +39 059 697705
sales.eu@elcam-medical.com

**China
Elcam China Office**

Room 3309, Jinhai Mansion
No. 89 of Anshun Road,
Shanghai, PRC. Zip: 200051
Tel: +86 21 62829967
Fax: +86 21 62829967
johnliu02@126.com



 Elcam Devices

 Elcam Stopcocks

 Elcam Accessories