



Precision & reliability at Every Point

Elcam Medical, a global leader in medical device manufacturing brings decades of OEM expertise to the bioprocessing suite.

Our integrated suite of components ensures integrity, safety, and precision throughout the fluid path.

Sensing Bio-P-Sensor™

Key features

Highly automated production line, 100%
leak testing in-production & Pre-calibration
Maximum pressure of 4 bar (60 psi)

Flow connections: ISO 80369-7 Standard
Luer, 1/8" and 3/16" barbs

Real-time
accuracy

High-precision pressure
monitoring to ensure process
stability and reliability



Customizable Flow Solutions

Multiple configurations
designed to streamline
your complex single-use
assemblies



Routing Bio-Manifolds

Key features

Integrated design: saving the need for additional
Luer to Barb adaptors

Customizable: optional color-coded handles
according to request

Flow connections: Elcam's 2-5 ports with luer
connection, 5 ports manifold with 1/8" or 3/16" barbs.

Closed, safe access

A closed-system design for protecting your process from contamination



Access Bio-NIP® NeedleFree Connector

Key features

Normally closed Luer activated valve, transparent fluid path, dedicated non-membrane activating cover cap for ensuring no air or microorganisms ingress

Flow connections: Elcam's off-the shelf NIP connectors; 1/8" or 3/16" barbs

Safety & Convenience Bio-Safe2 Rotator™

Key features

Featuring a unique flexible material (polyethylene) and spinning lock feature that provide ease of use and an extra layer of passive protection Ideal for cases where stable luer connection has to be made with barbed connection.

Flow connections: S2R with 1/8" and 3/16" barbs.

Connections you can trust

Minimizes leakages, prevents tube kinking and accidental disconnection with the unique rotating lock mechanism.



Precision Flow Control

Leveraging Elcam's global leadership to provide high-quality flow control



Flow Control Bio-Flow Valve World-Class Stopcocks


Key features

Efficiency: eliminating the need for additional Luer to Barb adaptors. High Variety – flow rates, handle types and colors according to request

Flow connections: Elcam's off-the-shelf standard Luer stopcocks variety, stopcocks with 1/8" and 3/16" barbs



Our integrated suite of components.

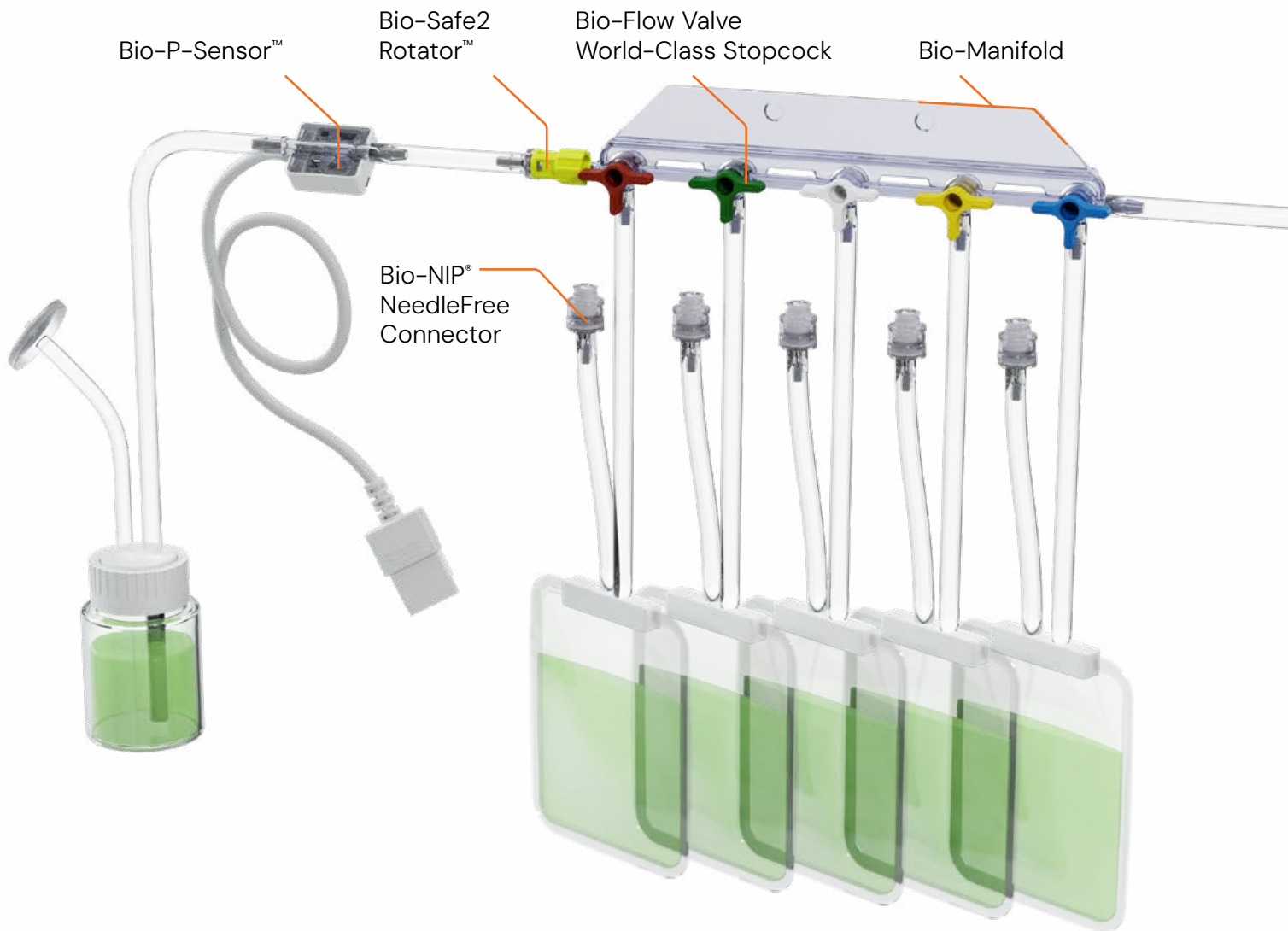

Samples available. Soon to be commercially available

elcam-medical.com

Head Office / Manufacturer Elcam Medical

BarAm 1386000, Israel
Tel: 972 4 698 8120/1
sales@elcam.co.il

USA Elcam Medical Inc.

1 University Plaza, Suite 511
Hackensack, NJ 07601
Tel: +1 201 4571120
info@elcam-medical.com

Europe Elcam Medical Europe

Via Emilia Romagna 15, 41012
Carpi (Modena) Italy
Tel: +39 059 697710
sales.eu@elcam-medical.com

China Elcam China Office

Room 3309, Jinhai Mansion
No. 89 of Anshun Road,
Shanghai, PRC. Zip: 200051
Tel: +86 21 62829967
johnliu02@126.com



Status: under development

**YOUR SOURCE OF CONSISTENT
QUALITY & CARE**