


Sense-IT™ (DIPT)

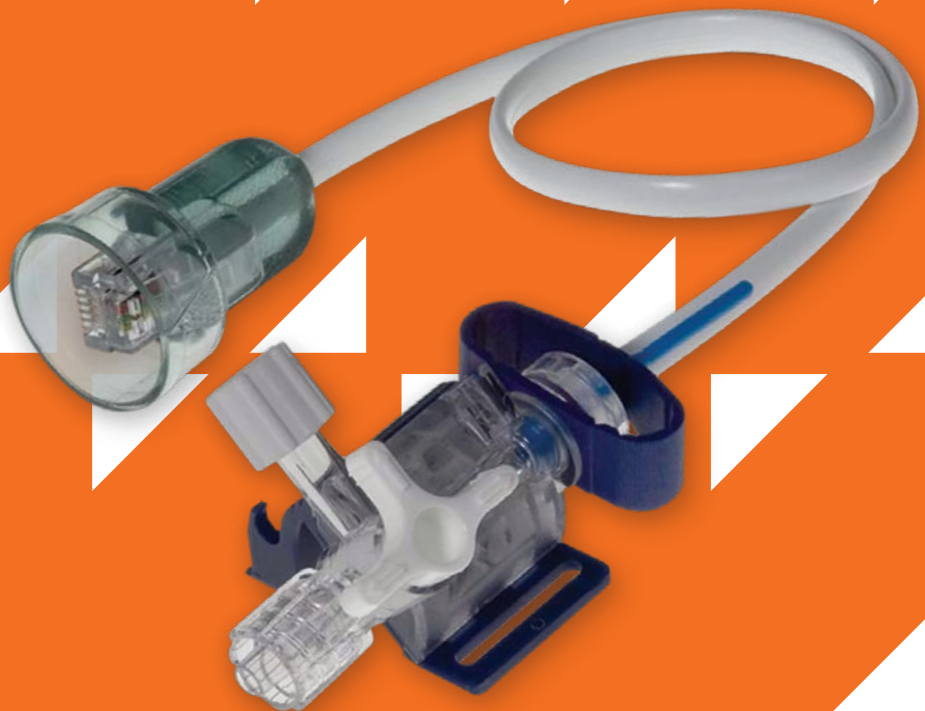

Elcam Disposable Integrated Pressure Transducer is indicated for pressure sensing during clinically invasive procedures.

Elcam Medical's Sense-IT™ is the first and only fully integrated pressure transducer. This integration provides you with the highest degree of reliability and robustness, along with the space saving and ease of use so much needed in Intensive Care Units and Operating Rooms.

Elcam offers a wide variety of configuration options for the different parts of the Sense-IT™

Benefits

Reliability

- Robust structural integrity
- Air-ingress safety
- Assured electrical and leakage conditions by 100% in-production tests
- State of the art assembly process
- Highly automated production line, fully controlled by on-line testing

Ease of use

- Full visibility of the entire flow path
- Easy and safe priming and debubbling
- Single hand operation of flush device
- Dual-mode flush device
- Compact size

Accuracy

- Very low dynamic response


Available configurations

Stopcock options: 3&4 way, handle shape and color, optional cap on female port (several venting and color options), tube assembly or tube fitment options.

3 or 30 ml/hr flush device: 30 ml version is designed especially for pediatric use, employing a pump, while the 3 ml/hr device is for more general use (gravity feed).

Connector type, Components color (base, elastomer), Printing of customer's logo: Available in a wide variation of options.

Available connectors:

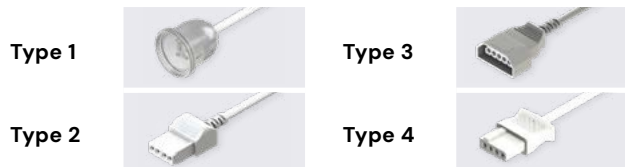

Applications

- Continuous arterial blood pressure monitoring
- Central venous pressure monitoring
- Pulmonary artery monitoring
- Left atrial monitoring when used with an air eliminator
- Cardiac catheterization

Sense-IT™ (DIPT) specification

Sensing unit properties	
Pressure Ranges	-30 to +300 mmHg (-0.58 to +5.8psi)
Overpressure	Min. 6465 mmHg (125 psi)
Operating Product Life	168 hours
Shelf Life	5 years
Performance Parameters ¹	
Frequency Response	Min. 1200 Hz
Offset Drift ²	Max. 1mmHg/ 8hrs
Mechanical specifications	
Low flow test (3ml/hr variant)	1.9 to 5.0 ml/hr@276mmHg
Low flow test (30ml/hr variant)	20 to 40 ml/hr@276mmHg
High flow test	Min. 125 ml/min@276mmHg
Multiple engagements	Continuity is kept following 30 cable engagements

¹ Assuming 6 VDC and 22°C and after 5 sec. warm-up

² Normalized offset/bridge voltage — 8 hours after 20 sec. warm-up

Regulatory status:

- 510(k) cleared
- Biocompatible according to ISO 10993-1
- Compliance with the following international standards:
ANSI/AAMI BP22-1994 IEC 60601-1, IEC 60601-2-34

Intended use:

For single use only.
Maximum use duration
is 96 hours.



www.elcam-medical.com

Head Office/Manufacturer Elcam Medical

BarAm 1386000,
Israel
Tel: 972 4 698 8120/1
Fax: 972 4 698 0777
sales@elcam.co.il

USA Elcam Medical Inc.

1 University Plaza, Ste. 511
Hackensack, NJ 07601
Tel: +1 201 4571120
Fax: +1 201 4571125
info@elcam-medical.com

Europe Elcam Medical Europe

Via Emilia Romagna 15,
41012 Carpi (Modena) Italy
Tel: +39 059 697710
Fax: +39 059 697705
sales.eu@elcam-medical.com

China Elcam China Office

Room 3309, Jinhai Mansion
No. 89 of Anshun Road,
Shanghai, PRC. Zip: 200051
Tel: +86 21 62829967
Fax: +86 21 62829967
johnliu02@126.com